

Flexineb[®]2


Testimonials

"We have used this product religiously leading up the biggest four star in the world."

I can say, for the first time in two years, and at the biggest, fastest, longest track of his career – Tate did not bleed."

Megan Kepferle

Head groom at Sinead Halpin Eventing

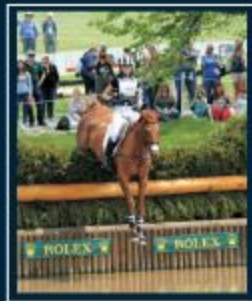

"The Flexineb is the best thing I have had as an aid to my horses! It is a crucially important piece of equipment and so easy to use."

Jessica Harrington


Leading Irish Dual Racehorse Trainer



"Regarding the Flexineb, I have used the adult mask in a few horses with pneumonia and I was impressed with how well it works and how silent it is. The horses tolerate it very well and it is much faster than what I had used before."

Steeve Giguere DVM, PhD, DACVIM

Professor and Hodgson Research Chair in Equine Studies
College of Veterinary Medicine
University of Georgia

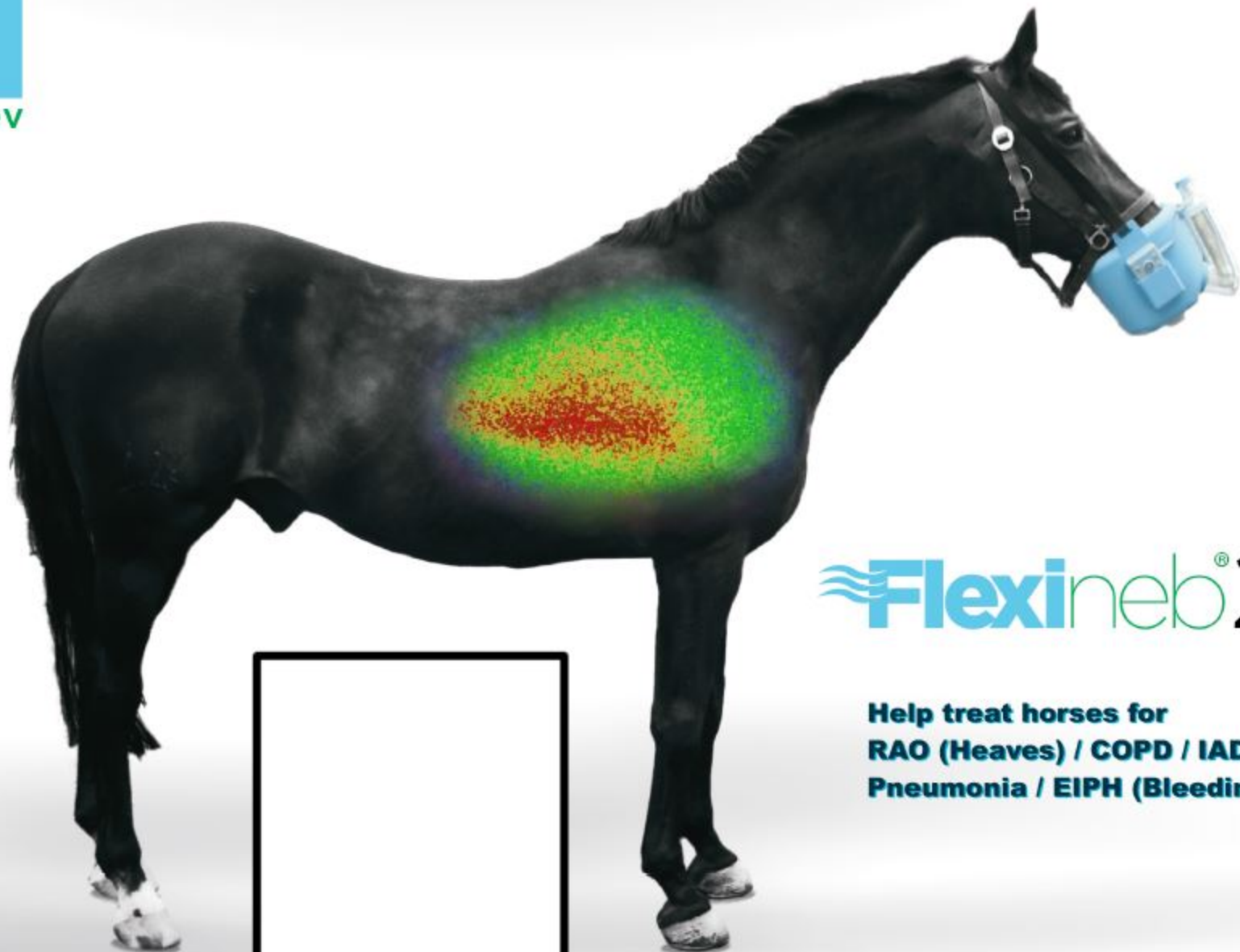

Flexineb[®]2

**Help treat horses for
RAO (Heaves) / COPD / IAD
Pneumonia / EIPH (Bleeding)**

What is nebulization?

Nebulization is a treatment technique in which tiny spherical droplets of aerosolised medication are generated and inhaled by the patient.

Treatment Method	Delivery Median
Inter-venous	Bloodstream
Oral - Tablet	Digestive system
Oral - Feed	Digestive system
Nebulizer	Direct to affected site, ie lungs



Why use a Nebulizer?

A nebulizer device enables the horse owner to easily administer care to their horses airways and lungs. The device incorporates a mask that fits on the nose and a technological means of generating a fine mist. Old style technologies are based on compressor technologies which plug into the mains supply. Newer technologies are ultrasonic devices, while the most up-to-date devices (e.g. Flexineb) incorporate a vibrating piezo component. Each of these technologies can generate the fine mist which ensures the mist goes directly into the airways and ultimately the lungs.

Nebulizers:

- Allow the treatment of breathing problems diagnosed by professional veterinarians.
- Allow easy conditioning of airways to help prepare for competition.

Applications:

- Deliver antibiotics for infections, corticosteroids for inflammation, bronchodilators to open airways.
- Administer saline and/or natural solutions to condition airways (e.g. remove mucus).
- Clear dirty noses.
- Preparation for sales.
- Manage dry coughs.
- Clear dust allergies.
- Treat pneumonia.

Please refer to testimonials on www.nortev.com to see how Flexineb users apply this product.

Why Flexineb?

Extremely lightweight (only 1kg) for comfort of horse and ease of use for caregiver

Fully integrated blue rubber mask, lithium ion rechargeable battery powered box, and medication holding cup - ALL included in system when shipped

Revolutionary blue flexible mask ensures maximum tolerance by the horse

Fastest delivery times on the market combined with excellent portability and silent operation

Designed and developed in Ireland with the help of Irish and International veterinarians, thoroughbred trainers, and horse owners.



Picture courtesy of Hagyard Equine Medical Institute



Flexineb[®]2


Our new Soft Mount Med Cup with refined membrane hole size recently tested at Melbourne Scientific has achieved the following results:

VMD Dv (50) 3.105 microns
 % <5 microns = 70.57%
 % 2-5 microns = 39.02%
 % 1 microns = 23.75%

- 1 Upper Airway (naso-pharynx) >10 μ
- 2 Central Airway (trachea) 10-5μ
- 3 Lower Airways (lungs) 5-2μ
- 4 Alveoli 1μ

NEW: Flexineb[®]2

- ≡ Longer Battery Life
- ≡ New Input Connector (moisture Protection)
- ≡ New Easier to Use Output Cable
- ≡ New AUTOMATIC TUNING (AMT)
Automatically adjusts vibrations to improve medicine cup function

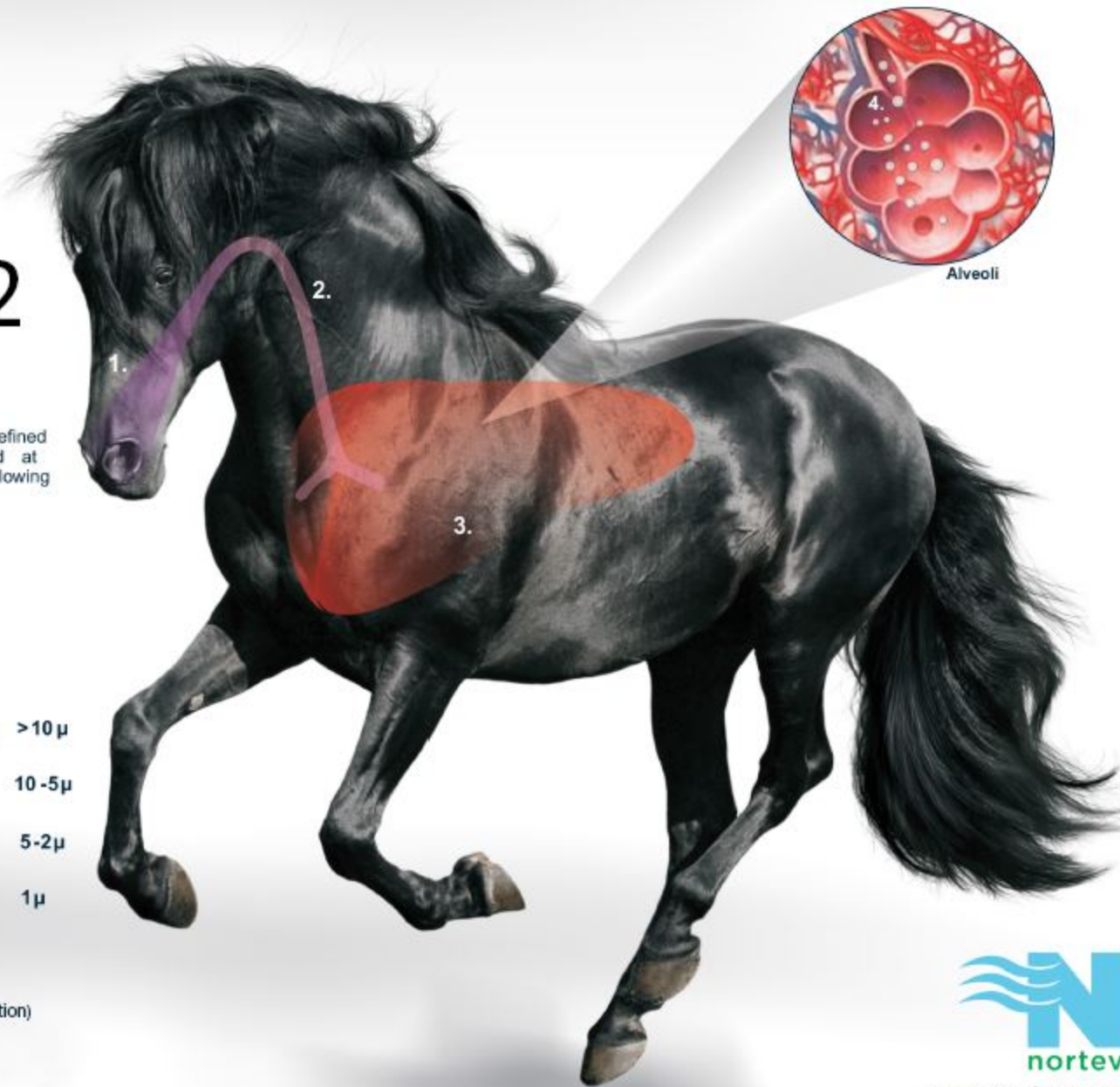

www.nortev.com